



Ordering Information

ASAHI Masters PARKWAY SOFT

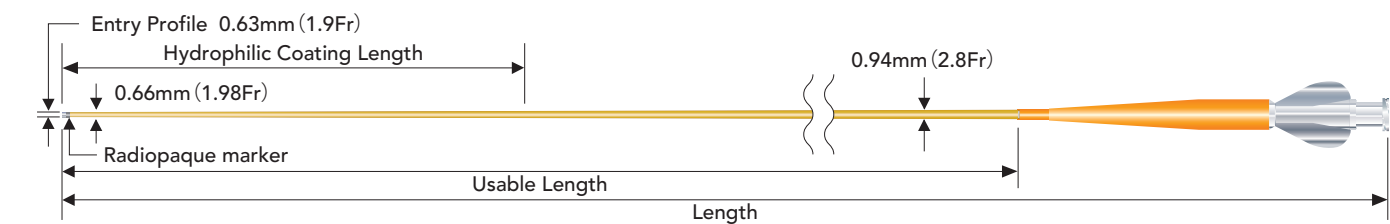


Product	Catalog No.	Usable Length	Length	O.D. of Distal shaft / O.D. of Proximal shaft	Shaft I.D.	Max. Guide Wire	Hydrophilic Coating Length	Tip Shape	Max. Pressure
ASAHI Masters PARKWAY SOFT	WMST105-18PWSF	105cm	112cm	0.66mm / 1.98Fr 0.94mm / 2.8Fr	0.55mm / 0.022inch	0.46mm / 0.018inch	65cm	Straight	6895kPa / 1000psi
	WMST125-18PWSF	125cm	132cm	0.66mm / 1.98Fr 0.94mm / 2.8Fr	0.55mm / 0.022inch	0.46mm / 0.018inch	65cm	Straight	6895kPa / 1000psi
	WMST150-18PWSF	150cm	157cm	0.66mm / 1.98Fr 0.94mm / 2.8Fr	0.55mm / 0.022inch	0.46mm / 0.018inch	100cm	Straight	6895kPa / 1000psi
	WMST45A-18PWSSF	105cm	112cm	0.66mm / 1.98Fr 0.94mm / 2.8Fr	0.55mm / 0.022inch	0.46mm / 0.018inch	65cm	Angled	6895kPa / 1000psi
	WMST45A-18PWSF	125cm	132cm	0.66mm / 1.98Fr 0.94mm / 2.8Fr	0.55mm / 0.022inch	0.46mm / 0.018inch	65cm	Angled	6895kPa / 1000psi

ASAHI Masters PARKWAY HF

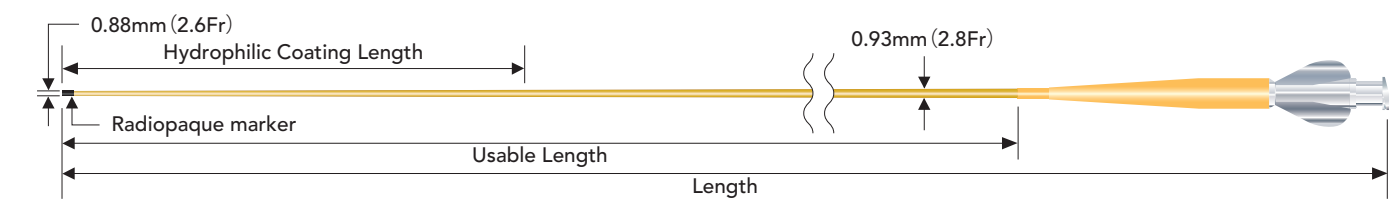

Product	Catalog No.	Usable Length	Length	O.D. of Distal shaft / O.D. of Proximal shaft	Shaft I.D.	Max. Guide Wire	Hydrophilic Coating Length	Tip Shape	Max. Pressure
ASAHI Masters PARKWAY HF	WMST105-27HF	105cm	112cm	0.88mm / 2.6Fr 0.93mm / 2.8Fr	0.69mm / 0.027inch	0.53mm / 0.021inch	65cm	Straight	6895kPa / 1000psi
	WMST125-27HF	125cm	132cm	0.88mm / 2.6Fr 0.93mm / 2.8Fr	0.69mm / 0.027inch	0.53mm / 0.021inch	65cm	Straight	6895kPa / 1000psi
	WMST150-27HF	150cm	157cm	0.88mm / 2.6Fr 0.93mm / 2.8Fr	0.69mm / 0.027inch	0.53mm / 0.021inch	100cm	Straight	6895kPa / 1000psi

ASAHI Masters PARKWAY HF KIT

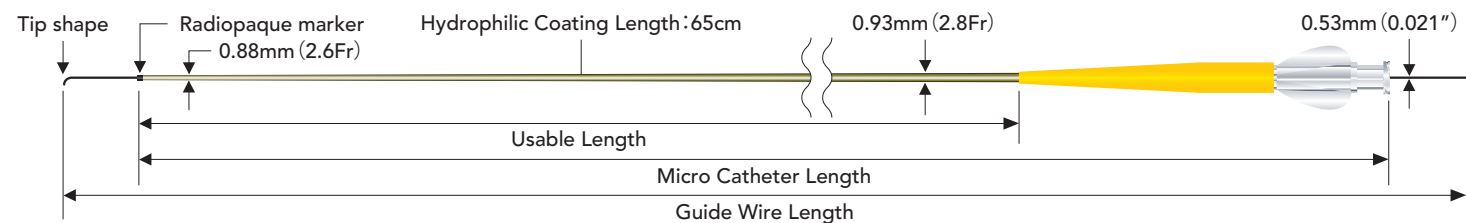

Product	Catalog No.	Micro Catheter Usable Length	Micro Catheter Length	O.D. of Distal shaft / O.D. of Proximal shaft	Shaft I.D.	Guide Wire Shaft O.D.	Guide Wire Length	Tip shape of GuideWire	Max. Pressure
ASAHI Masters PARKWAY HF KIT	WMST105-27HFK	105cm	112cm	0.88mm / 2.6Fr 0.93mm / 2.8Fr	0.69mm / 0.027inch	0.53mm / 0.021inch	140cm	Angled 80°	6895kPa / 1000psi
	WMST125-27HFK	125cm	132cm	0.88mm / 2.6Fr 0.93mm / 2.8Fr	0.69mm / 0.027inch	0.53mm / 0.021inch	160cm	Angled 80°	6895kPa / 1000psi
	WMST105-27HFKM	105cm	112cm	0.88mm / 2.6Fr 0.93mm / 2.8Fr	0.69mm / 0.027inch	0.53mm / 0.021inch	140cm	Multi Curve	6895kPa / 1000psi
	WMST125-27HFKM	125cm	132cm	0.88mm / 2.6Fr 0.93mm / 2.8Fr	0.69mm / 0.027inch	0.53mm / 0.021inch	160cm	Multi Curve	6895kPa / 1000psi



ASAHI and MastersPARKWAY are trademarks or registered trademarks of ASAHI INTECC CO., LTD. in Japan and other countries.

Manufactured by:



ASAHI INTECC CO.,LTD.
<http://www.asahi-intecc.com>

Distributed by:

ASAHI Masters PARKWAY SOFT
ASAHI Masters PARKWAY HF
ASAHI Masters PARKWAY HF KIT

Microcatheter / Guide wire Pre-loaded Microcatheter




ASAHI MastersPARKWAY series/ver.1/131225

ASAHI Masters PARKWAY SOFT

Microcatheter

- Radius tip and small profile for distal accessibility.
- Large ID for high volume injection and smooth delivery of embolic materials.
- Tungsten braiding for lumen integrity and visibility.
- Pre-shaped tip.

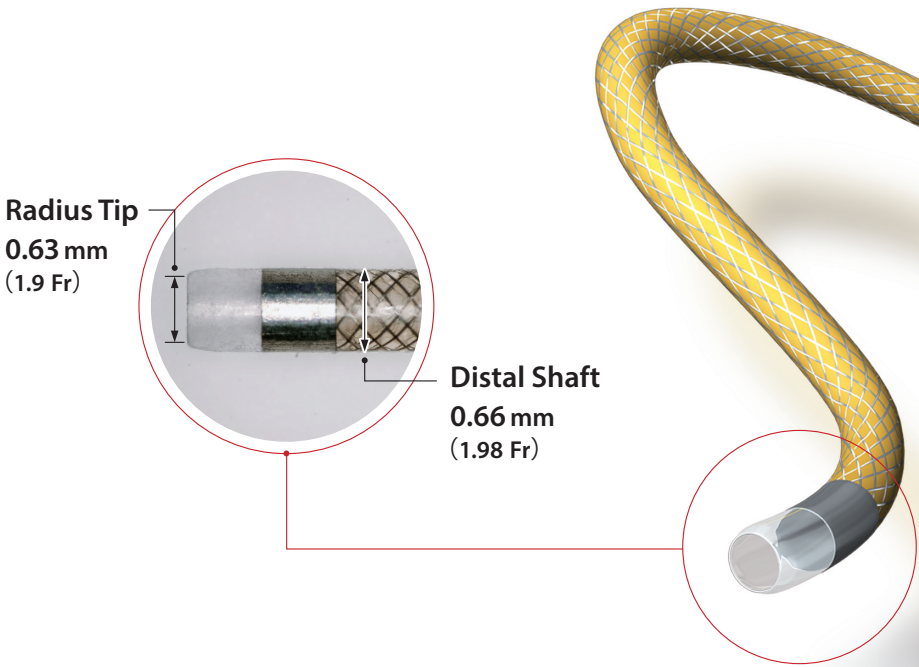

Flow Rate

Product	Usable Length	Contrast	Iodine Content	Flow Rate	Dead Volume
ASAHI Masters PARKWAY SOFT	105cm	Iopamidol	300 (mg/ml)	2.6 (ml/sec) / 1000 psi	0.43ml
			370 (mg/ml)	1.5 (ml/sec) / 1000 psi	
	125cm		300 (mg/ml)	2.0 (ml/sec) / 1000 psi	0.48ml
			370 (mg/ml)	1.0 (ml/sec) / 1000 psi	
	150cm		300 (mg/ml)	1.6 (ml/sec) / 1000 psi	0.54ml
			370 (mg/ml)	0.8 (ml/sec) / 1000 psi	

The above data was obtained by company standardized test, which may differ from industry standardized tests.
The above data does not prove that all devices have exactly the same performance with the samples used for these tests.


ASAHI Masters PARKWAY HF

Microcatheter

ASAHI Masters PARKWAY HF KIT

Guide wire Pre-loaded Microcatheter

- Tungsten braiding for lumen integrity and visibility.
- High volume of contrast and therapeutic agents through 0.027 inch (0.69 mm) large inner-lumen.
- 0.021 inch (0.53 mm) guide wire pre-loaded for Masters PARKWAY HF kits.

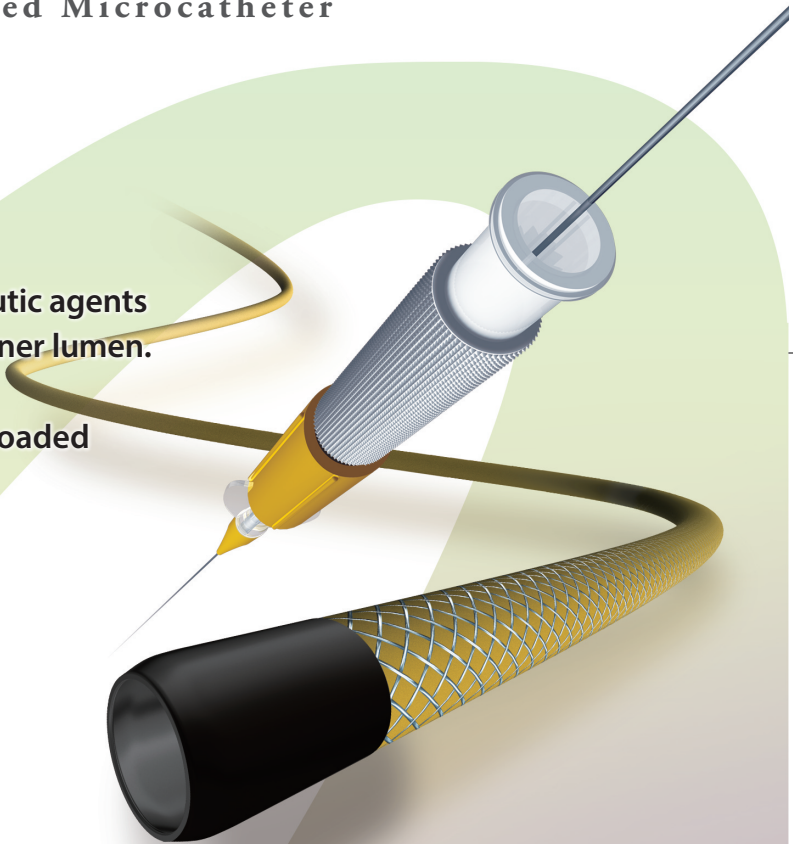

Flow Rate

Product	Usable Length	Contrast	Iodine Content	Flow Rate	Dead Volume
ASAHI Masters PARKWAY HF	105cm	Iopamidol	300 (mg/ml)	5.1 (ml/sec) / 1000psi	0.62ml
			370 (mg/ml)	2.7 (ml/sec) / 1000psi	
	125cm		300 (mg/ml)	4.4 (ml/sec) / 1000psi	0.69ml
			370 (mg/ml)	2.4 (ml/sec) / 1000psi	
	150cm		300 (mg/ml)	3.8 (ml/sec) / 1000psi	0.79ml
			370 (mg/ml)	2.0 (ml/sec) / 1000psi	

The above data was obtained by company standardized test, which may differ from industry standardized tests.
The above data does not prove that all devices have exactly the same performance with the samples used for these tests.