
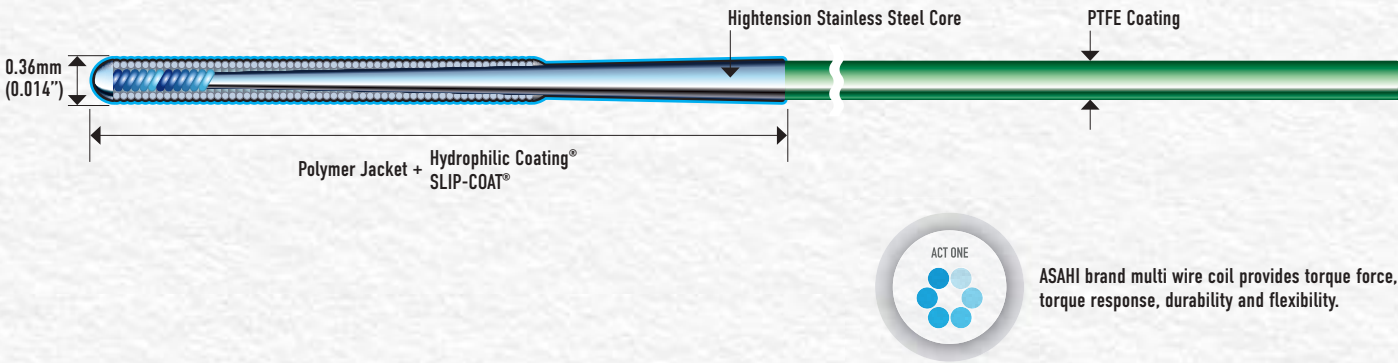


Structure & Ordering Information								
Product	Catalog No.	Diameter	Usable Length	Spring Coil Length	Coating Length	Radiopaque Length	Tip Shape	Clip Color
ASAHI Gladius MG	AP14R023P	0.36mm (0.014 inch)	190cm	8.5cm	SLIP-COAT® 41cm	3cm	Pre-shape	
	AP14R323P		300cm					



ASAHI Gladius MG



ASAHI Gladius MG

PCI GUIDE WIRE



Some materials/graphics were designed & produced by Vascular Perspective Ltd. www.vascularperspectives.com.

© 2022 ASAHI INTECC CO., LTD.
"ASAHI" and "Gladius" are trademarks or registered trademarks of ASAHI INTECC CO., LTD. in Japan and other countries. Coated with SLIP-COAT® coating

Manufactured by: Contact:

ASAHI INTECC EUROPE B.V.

ASAHI INTECC Europe B.V.
Strawinskylaan 527, 1077 XX Amsterdam, The Netherlands
Contact phone number: +31-20-794-0640
www.asahi-intecc.com



for Android for iOS
ASAHI product information
K22207_B_LC_2.0

YOUR DREAMS. WOVEN TOGETHER.


ASAHI Gladius MG

Broaden your approach to CTO cases



Characteristics


- Narrow loop and durability delivered through the modified distal core design
- High lubricity for smooth manipulation
- High torque response provided by ACT ONE technology.

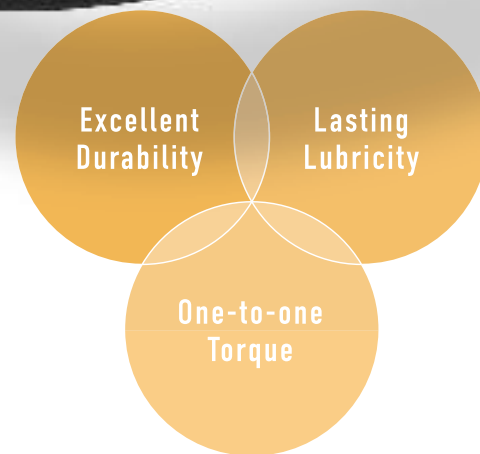


Excellent Durability

ASAHI Gladius MG maintains shape and guidewire performance even after crossing the lesion.


Guidewire prolapse examination after lesion crossing simulation test

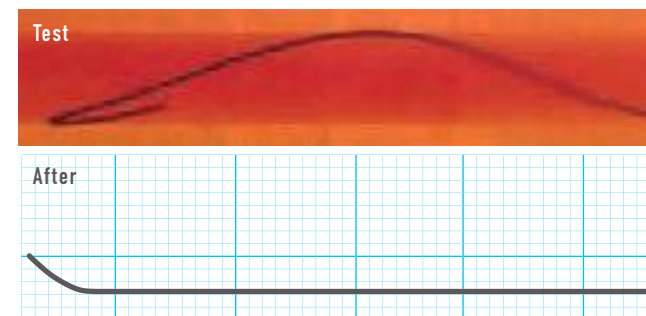

ASAHI Gladius MG

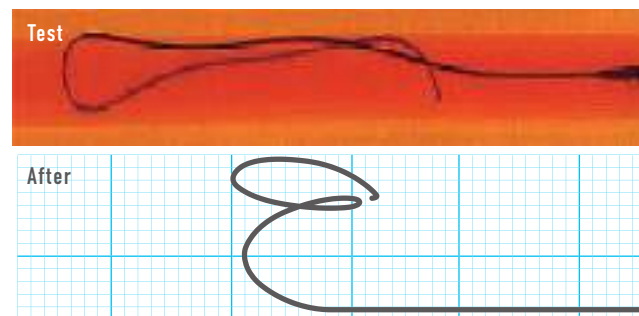


Competitor A

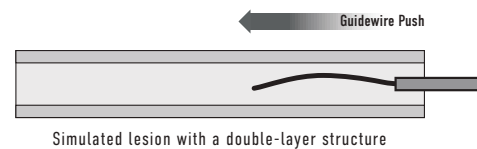


Torque Response Test

ASAHI Gladius MG with ACT ONE provides quick initial torque and superior torque response without whip motion.

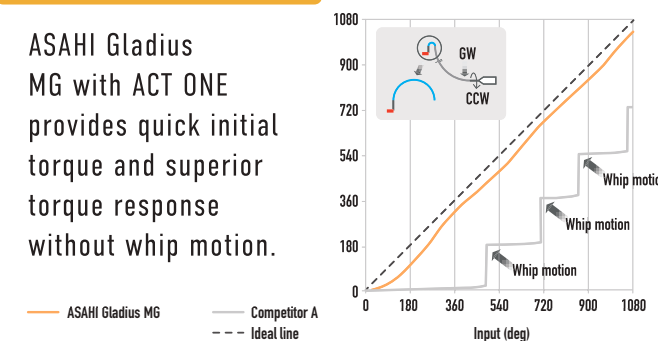



Lubricity Test

ASAHI Gladius MG has excellent lubricity that allows for smooth GW manipulation.

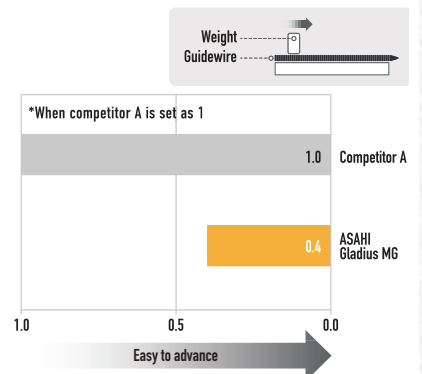

Flexibility Test

ASAHI Gladius MG has a modified distal core balance that minimizes a prolapse occurring in the tip.

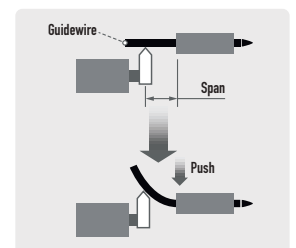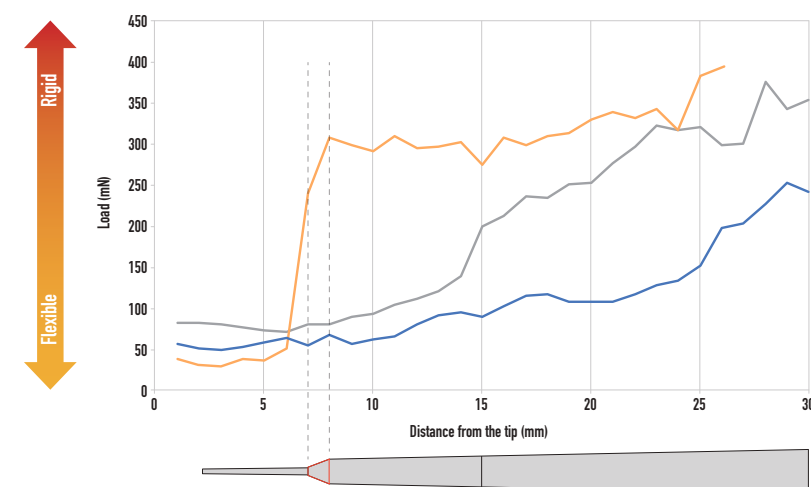

ASAHI Gladius EX
ASAHI Gladius MG
Competitor A